



Genotyping protocol

Iqce

IR00003909 / E258

(ICS internal reference)

This report has been prepared by: **David MOULAERT**
genotyping@igbmc.fr

This report has been validated by: **Sylvie Jacquot, PhD, Head of Genotyping Service**
33 (0)3 88 65 57 44
genotyping@igbmc.fr

The first version of this report was generated the: 23 Feb 2016

For any question, please contact:

Institut Clinique de la Souris - ICS - Mouse Clinical Institute
1 rue Laurent Fries, BP 10142
67404 Illkirch Cedex, France
Email: genotyping@igbmc.fr
Web site: <http://www-mci.u-strasbg.fr/>

TABLE OF CONTENTS

Table of contents2

1. Genotyping protocol and data2

 1.1. Genotyping strategy2

 1.2. PCR protocol4

2. Cre and Flp genotyping method5

1. Genotyping protocol and data


This section describes the condition used at the Mouse Clinical Institute (ICS) to genotype your Iqce Constitutive Knockout / Conditional Knockout (KO-cKO x Cre) project.

1.1. Genotyping strategy

The map below describes the position of the primers used for genotyping for each possible allele.



KO-cKO pc _ Genotyping strategy

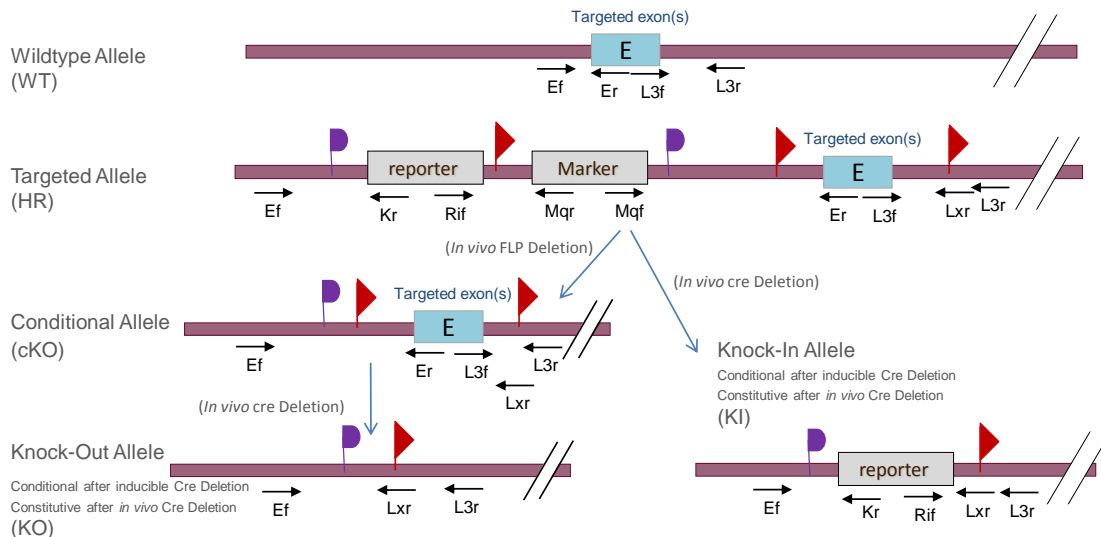

Sequence of primers used for genotyping:

Position	Primers	Sequence
Ef	6694	GAAAGGCAGGTCCTCTCTATAGC
Ef ²	6693	CCTAGAAGCTGCAGGACCTTGTTGC
Er	6698	CTGTGTGTCCACAGCAGATGCACAAAGAG
Kr	3277	CTCCTACATAGTTGGCAGTGTGGG
L3f	6696	CACAAACCTCCTCCACATCAC
L3f ²	6695	CACATCACCAAGTAAGTATGGCTGGAC
L3r	6697	GCATGACAACACAGCTTAAGGGCTTC
Lxr	3255	ACTGATGGCGAGCTCAGACCATAAC
Ri1f	5966	GCACATGGCTGAATATCGACGGT

²: for a selected position, a second primer was designed

PCR fragments expected size (bp):

Region analyzed	Primers used	Position on the primer (see the map above)	Targeted allele (HR)	conditional allele (cKO)	KI allele	WildType allele
5' part of the selection marker	6694-3277	Ef / Kr	327	---	327	---
Presence of the distal loxP	6696-6697	L3f / L3r	382	382	---	395
Distal loxP specific PCR	6695-3255	L3f ² / Lxr	251	251	---	---
Excision of the selection marker	6693-6698	Ef ² / Er	7283*	365	---	170
Cre total excision	5966-3255	Ri1f / Lxr	3133*	---	471	---

*: this PCR product will not be observed using our PCR genotyping conditions (see description below)

** : this PCR is only verified if mice are generated

---: no Amplicon should be obtained

1.2. PCR protocol

This section describes the composition of the mix and cycling conditions used for genotyping.

Reagents:	Volume:
- FastStart PCR Master (Roche)	7.5µl
- DNA (50ng/µl)	1.5µl
- 5' primer (100 µM)	0.06µl
- 3' primer (100 µM)	0.06µl
- Sterile H ₂ O	up to 15 µl

Cycling conditions:

Temp	Time	#Cycles
95°C	4min	1
94°C	30s	34
62°C	30s	
72°C	1min	
72°C	7min	1
20°C	5min	1

NB: These PCR conditions have been optimized for high-throughput genotyping. Adaptation to small-scale may be required.

2. Cre and Flp genotyping method

You will find the genotyping protocol in the publication:

[Highly-efficient, fluorescent, locus directed cre and FlpO deleter mice on a pure C57BL/6N genetic background.](#)

Birling MC, Dierich A, Jacquot S, Hérault Y, Pavlovic G.
Genesis. 2012 Jun;50(6):482-9. doi: 10.1002/dvg.20826. Epub 2012 Mar 20.