



Genotyping protocol
Lcp2tm1a(KOMP)Wtsi
EPD0396_3_B12

This report has been prepared by: **Valérie Rousseau**
genotyping@igbmc.fr

This report has been validated by: **Sylvie Jacquot, PhD, Head of Genotyping Service**
33 (0)3 88 65 57 44
genotyping@igbmc.fr

The first version of this report was generated the: 05 Dec 2014

For any question, please contact:

Institut Clinique de la Souris - ICS - Mouse Clinical Institute
1 rue Laurent Fries, BP 10142
67404 Illkirch Cedex, France
Email: genotyping@igbmc.fr
Web site: <http://www-mci.u-strasbg.fr/>

TABLE OF CONTENTS

Table of contents2

1. Genotyping protocol and data2

 1.1. Genotyping strategy2

 1.2. PCR protocol4

 1.3. Picture of genotyping with various alleles5


2. Cre and Flp genotyping method6

1. Genotyping protocol and data

This section describes the condition used at the Mouse Clinical Institute (ICS) to genotype your **Lcp2** Constitutive Knockout / Conditional Knockout (KO-cKO x Cre) project.

1.1. Genotyping strategy

The map below describes the position of the primers used for genotyping for each possible allele.

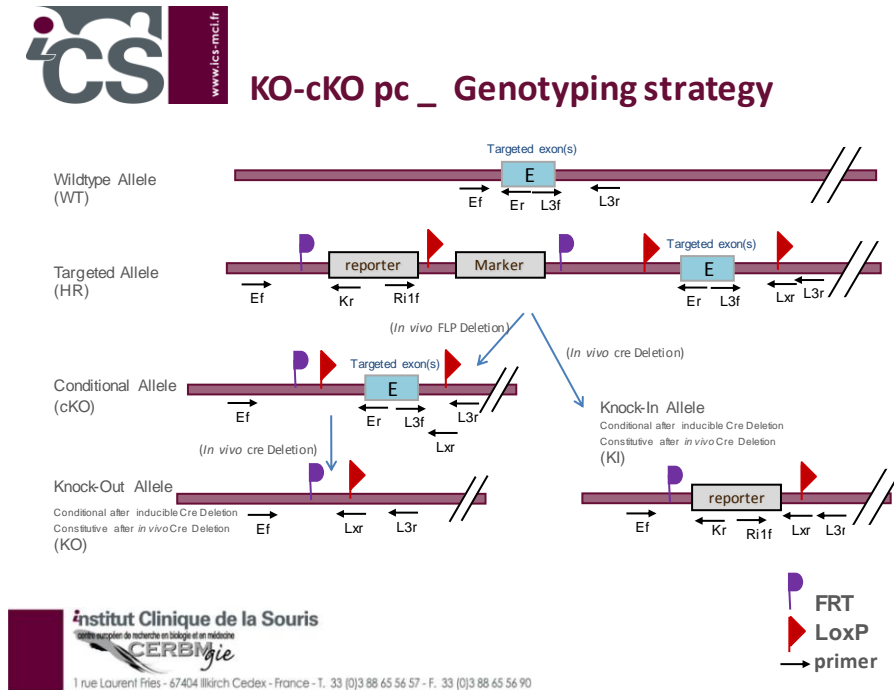

Sequence of primers used for genotyping:

Position	Primers	Sequence
Ef	7530	CACCTTTCAGAGTCCAATGGACTATGC
Er	7533	CTAATAAGACAAGGGACAATCTG
Kr	3277	CTCCTACATAGTTGGCAGTGTGGG
L3f	7531	GGTGGACGATGCTGCTGACTATGA
L3r	7532	GCACATGGTGTCTGGCTCATAATAGGA
Lxr	3254	TTATCATTAAATTGCGTTGCGCCATC
Lxr2	3255	ACTGATGGCGAGCTCAGACCATAAC
Ri1f	5966	GCACATGGCTGAATATCGACGGT

²: for a selected position, a second primer was designed

PCR fragments expected size (bp):

Region analyzed	Primers used	Position on the primer (see the map above)	Targeted allele (HR)	conditional allele (cKO)	KI allele	WildType allele
5' part of the selection marker	7530-3277	Ef / Kr	395	---	---	---
Presence of the distal loxP	7531-7532	L3f / L3r	475	475	---	484
Distal loxP specific PCR	7531-3254	L3f / Lxr	305	305	---	---
Excision of the selection marker	7530-7533	Ef / Er	7370*	466	---	287
Cre total excision	5966-3255	Ri1f / Lxr2	---*	---	471	---

*: this PCR product will not be observed using our PCR genotyping conditions (see description below)

---: no Amplicon should be obtained

1.2. PCR protocol

This section describes the composition of the mix and cycling conditions used for genotyping.

Reagents:	Volume:
- FastStart PCR Master (Roche)	7.5µl
- DNA (50ng/µl)	1.5µl
- 5' primer (100 µM)	0.06µl
- 3' primer (100 µM)	0.06µl
- Sterile H ₂ O	up to 15 µl

Cycling conditions:


Temp	Time	#Cycles
95°C	4min	1
94°C	30s	34
62°C	30s	
72°C	1min	
72°C	7min	1
20°C	5min	1

NB: These PCR conditions have been optimized for high-throughput genotyping. Adaptation to small-scale may be required.

1.3. Picture of genotyping with various alleles

Analysis of PCR products pattern was done by gel electrophoresis 2% agarose (SB buffer).

Representative genotyping picture

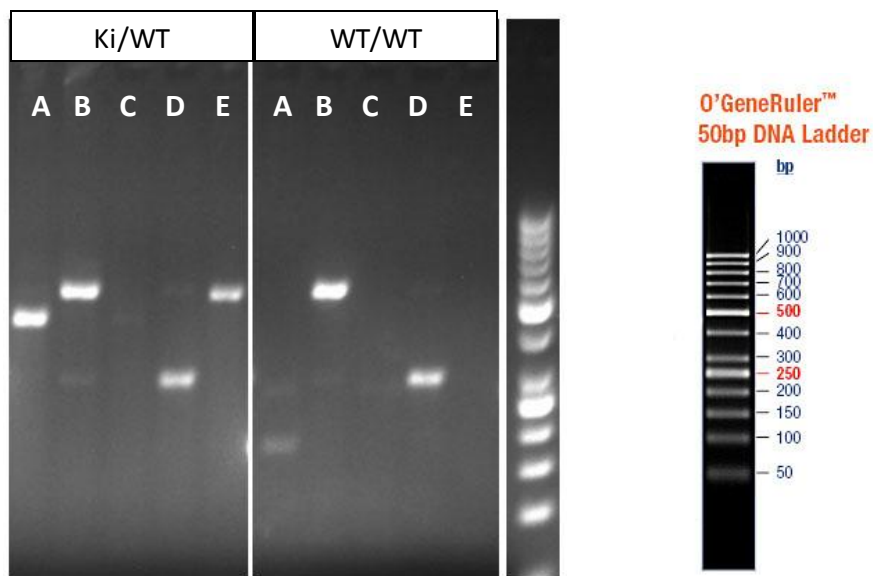

- A 5' part of the selection marker
- B Presence of the distal loxP
- C Distal loxP specific PCR
- D Excision of the selection marker
- E Cre total excision

2. Cre and Flp genotyping method

You will find the genotyping protocol in the publication:

[Highly-efficient, fluorescent, locus directed cre and FlpO deleter mice on a pure C57BL/6N genetic background.](#)

Birling MC, Dierich A, Jacquot S, Hérault Y, Pavlovic G.

Genesis. 2012 Jun;50(6):482-9. doi: 10.1002/dvg.20826. Epub 2012 Mar 20.