



Genotyping protocol
Stat1^{tm1a}(EUCOMM)Wtsi
EPD0496_3_D07

This report has been prepared by: **Sylvie Jacquot**
genotyping@igbmc.fr

This report has been validated by: **Sylvie Jacquot, PhD, Head of Genotyping Service**
33 (0)3 88 65 57 44
genotyping@igbmc.fr

The first version of this report was generated the: 7 nov 2014

For any question, please contact:

Institut Clinique de la Souris - ICS - Mouse Clinical Institute
1 rue Laurent Fries, BP 10142
67404 Illkirch Cedex, France
Email: genotyping@igbmc.fr
Web site: <http://www-mci.u-strasbg.fr/>

TABLE OF CONTENTS

Table of contents2

1. Genotyping protocol and data2

 1.1. Genotyping strategy2

 1.2. PCR protocol4


2. Cre and Flp genotyping method5

1. Genotyping protocol and data

This section describes the condition used at the Mouse Clinical Institute (ICS) to genotype your **Stat1** Constitutive Knockout / Conditional Knockout (KO-cKO x Cre) project.

1.1. Genotyping strategy

The map below describes the position of the primers used for genotyping for each possible allele.

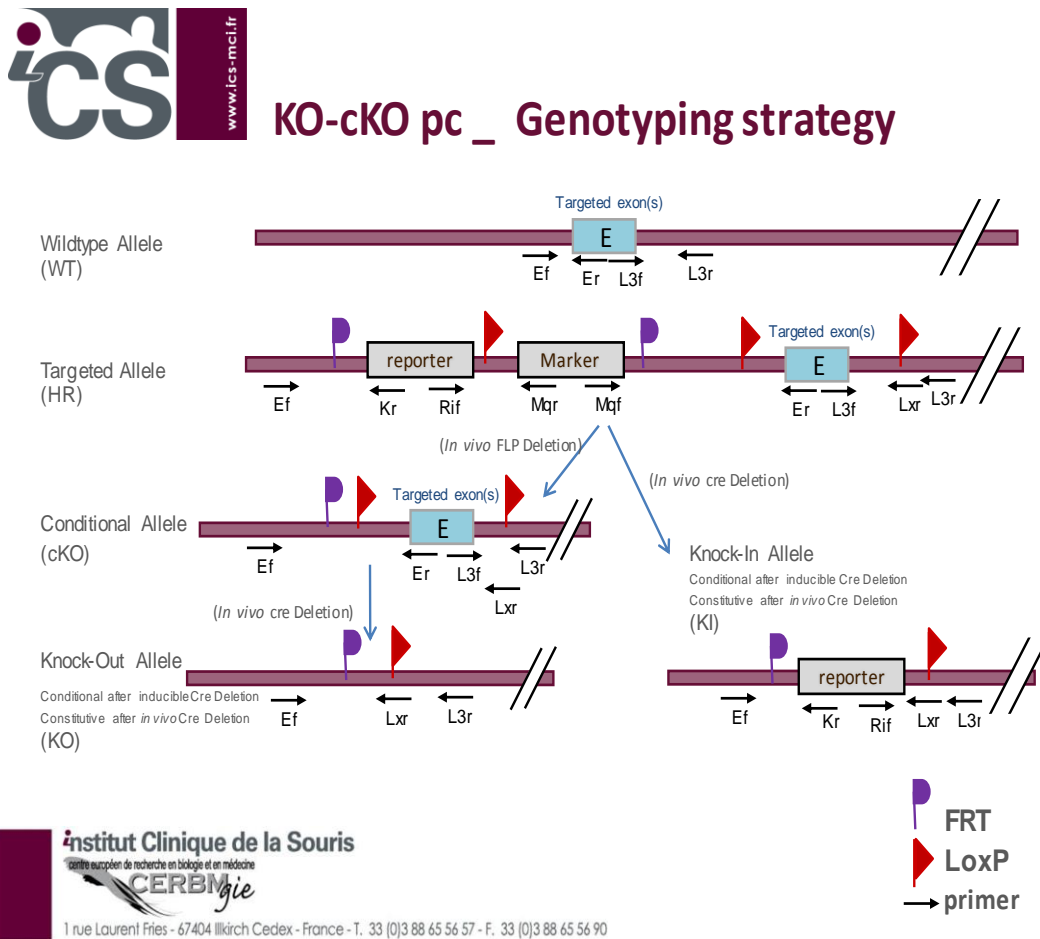

Sequence of primers used for genotyping:

Position	Primers	Sequence
Ef	8002	GGATATCCCATTCATCTCCTGAGTGC
Ef ²	8003	GATCAGTCTAGATCATCCGCCGTC
Er	8007	CCAGCATCCTGACCTTATGGAACACTAC
Kr	3209	CCAACAGCTTCCCCACAACGG
L3f	8005	CCTCTCAGAACAGAGGTAAGAGTCC
L3f ²	8004	GGCTGCTTTTTCCCCTGTGTACG
L3r	8006	CAGCGACAATAAAACACAATGCTGCAGC
Lxr	3255	ACTGATGGCGAGCTCAGACCATAAC
Rif	5966	GCACATGGCTGAATATCGACGGT

²: for a selected position, a second primer was designed

PCR fragments expected size (bp):

Region analyzed	Primers used	Position on the primer (see the map above)	Targeted allele (HR)	conditional allele (cKO)	KI allele	WildType allele
5' part of the selection marker	8002-3209	Ef / Kr	346	---	---	---
Presence of the distal loxP	8005-8006	L3f / L3r	355	355	---	387
Distal loxP specific PCR	8004-3255	L3f ² / Lxr	217	217	---	---
Excision of the selection marker	8003-8007	Ef ² / Er	7235*	331	---	210
Cre total excision	5966-3255	Rif / Lxr	3119*	---	471**	---

*: this PCR product will not be observed using our PCR genotyping conditions (see description below)

** : this PCR is only verified if mice are generated

---: no Amplicon should be obtained

1.2. PCR protocol

This section describes the composition of the mix and cycling conditions used for genotyping.

Reagents:	Volume:
- FastStart PCR Master (Roche)	7.5µl
- DNA (50ng/µl)	1.5µl
- 5' primer (100 µM)	0.06µl
- 3' primer (100 µM)	0.06µl
- Sterile H ₂ O	up to 15 µl

Cycling conditions:

Temp	Time	#Cycles
95°C	4min	1
94°C	30s	34
62°C	30s	
72°C	1min	
72°C	7min	1
20°C	5min	1

NB: These PCR conditions have been optimized for high-throughput genotyping. Adaptation to small-scale may be required.

2. Cre and Flp genotyping method

You will find the genotyping protocol in the publication:

[Highly-efficient, fluorescent, locus directed cre and FlpO deleter mice on a pure C57BL/6N genetic background.](#)

Birling MC, Dierich A, Jacquot S, Hérault Y, Pavlovic G.
Genesis. 2012 Jun;50(6):482-9. doi: 10.1002/dvg.20826. Epub 2012 Mar 20.