



Genotyping protocol

Mtmr12

IR00004493 / P4493

(ICS internal reference)

This report has been prepared by: **David MOULAERT**
genotyping@igbmc.fr

This report has been validated by: **Sylvie Jacquot, PhD, Head of Genotyping Service**
33 (0)3 88 65 57 44
genotyping@igbmc.fr

The first version of this report was generated the: 25 Feb 2016

For any question, please contact:

Institut Clinique de la Souris - ICS - Mouse Clinical Institute
1 rue Laurent Fries, BP 10142
67404 Illkirch Cedex, France
Email: genotyping@igbmc.fr
Web site: <http://www-mci.u-strasbg.fr/>

TABLE OF CONTENTS

Table of contents2

1. Genotyping protocol and data2

 1.1. Genotyping strategy2

 1.2. PCR protocol4

2. Cre and Flp genotyping method5

1. Genotyping protocol and data


This section describes the condition used at the Mouse Clinical Institute (ICS) to genotype your **Mtmr12** Constitutive Knockout / Conditional Knockout (KO-cKO x Cre) project.

1.1. Genotyping strategy

The map below describes the position of the primers used for genotyping for each possible allele.



KO-cKO pc _ Genotyping strategy

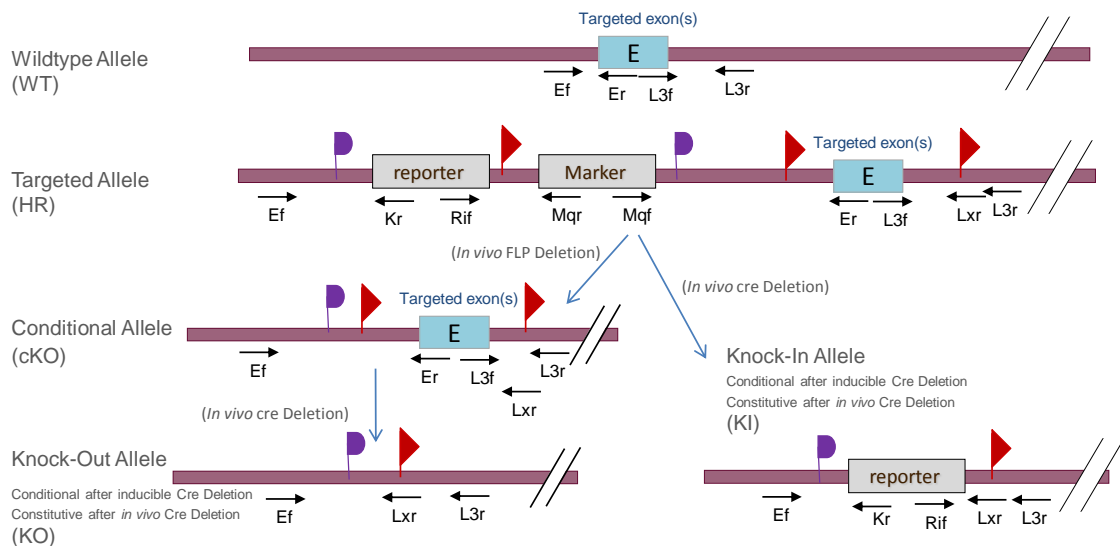

Sequence of primers used for genotyping:

Position	Primers	Sequence
Ef	6989	CTTACTGAGCCCAGGCTTACATGTCTG
Er	6994	GGTCTTTAGGCTTGGTAGCAAATAC
Er ²	6993	GAGGCTGGATTTCTCTTCCCACC
Kr	3278	GGGCAAGAACATAAAAGTGACCCTCC
L3f	6991	GAAGTCACTCTTCACTTGTTGCCAGG
L3f ²	6990	CGGCCATGTTATCAAGTCTCACATATAG
L3r	6992	CATGGGACGTTCTGCATTGGAGAC
Lxr	3255	ACTGATGGCGAGCTCAGACCATAAC
Ri1f	5966	GCACATGGCTGAATATCGACGGT

²: for a selected position, a second primer was designed

PCR fragments expected size (bp):

Region analyzed	Primers used	Position on the primer (see the map above)	Targeted allele (HR)	conditional allele (cKO)	KI allele	WildType allele
5' part of the selection marker	6989-3278	Ef / Kr	368	---	368	---
Presence of the distal loxP	6991-6992	L3f / L3r	343	343	---	344
Distal loxP specific PCR	6990-3255	L3f ² / Lxr	229	229	---	---
Excision of the selection marker	6989-6994	Ef / Er	7343*	439	---	335
Cre total excision	5966-3255	Ri1f / Lxr	3084*	---	471**	---

*: this PCR product will not be observed using our PCR genotyping conditions (see description below)

** : this PCR is only verified if mice are generated

---: no Amplicon should be obtained

1.2. PCR protocol

This section describes the composition of the mix and cycling conditions used for genotyping.

Reagents:	Volume:
- FastStart PCR Master (Roche)	7.5µl
- DNA (50ng/µl)	1.5µl
- 5' primer (100 µM)	0.06µl
- 3' primer (100 µM)	0.06µl
- Sterile H ₂ O	up to 15 µl

Cycling conditions:

Temp	Time	#Cycles
95°C	4min	1
94°C	30s	34
62°C	30s	
72°C	1min	
72°C	7min	1
20°C	5min	1

NB: These PCR conditions have been optimized for high-throughput genotyping. Adaptation to small-scale may be required.

2. Cre and Flp genotyping method

You will find the genotyping protocol in the publication:

[Highly-efficient, fluorescent, locus directed cre and FlpO deleter mice on a pure C57BL/6N genetic background.](#)

Birling MC, Dierich A, Jacquot S, Hérault Y, Pavlovic G.
Genesis. 2012 Jun;50(6):482-9. doi: 10.1002/dvg.20826. Epub 2012 Mar 20.