

EMMA ID: 13626
Gene: *Slc35f5*
Common name: *Slc35f5-em1_9*
Allele: *Slc35f5^{em1(IMPC)Hmgu}*

Allele Information

For more information on production, guides and mutation, search for gene/project, go to project summary, go to production plan, go to production outcome and "more details"

<https://www.gentar.org>

IMPC mouse phenotype data, search by the gene name
<http://www.mousephenotype.org/>

Genotyping Information

Genotyping by end-point PCR based on gel is composed of a genespecific short range PCR using primers on wild type allele and a mutant allele-specific short range PCR. The combined results show the genotype of the mice. For example: mutant positive, wild type positive = Heterozygous. In addition to the expected product, the mutant assay may also amplify the endogenous wild type sequence, which will appear as a larger band on an agarose gel. The presence of this extra band will depend on the size of the original deletion.

PCR primer pairs and expected size bands

Assay	Forward Primer	Reverse Primer	Expected Size Band (bp)
Wild type	Slc35f5_F	Slc35f5_wtR	501/1766
Mutant	Slc35f5_F	Slc35f5_R	317

Primer sequences

Primer Name	Sequence 5' --> 3'
Slc35f5_F	cctccaccaaatagttgtgttg
Slc35f5_R	cacaattgaaattgtggctccac
Slc35f5_wtR	ctcaacactccctcaacaactac

PCR setup (LongAmp® TaqDNA Polymerase)

Component	Volume (µl) 1x
DNA (~ 50-100 ng)	2-4
100% DMSO	0,4
PCR-Buffer (5x)	4
DNTP mix (10 mM)	0,5
Primer 1 (10 pmol/µl)	1
Primer 2 (10 pmol/µl)	1
Primer 3 (10 pmol/µl)	1
Taq Polymerase (2,5U/µl)	0,5
H2O*	7,6
Final volume	20


* The amount of H2O is adjusted with the number of primer.

Amplification conditions

PCR Settings	Temperature (°C)	Time	# of cycles
1 Denaturation (Melting)	94°C	3 min	1
2 Amplification (Melting, Annealing, Polym.)	94°C	30 sec	39
	65-55 (↓1°C/Cycle)	20 sec	
	65°C	1 min	
3 Polymerisation	65°C	10 min	1
4 Cooling	4°C	hold	1

Touch-Down cycling protocol: first 10 cycles anneal at 65°C, decreasing 1°C per cycle, next 30 cycles anneal at 55°C. These PCR conditions have been optimized for our methods and preparation kits. Adaptions may be required.

Gel Image

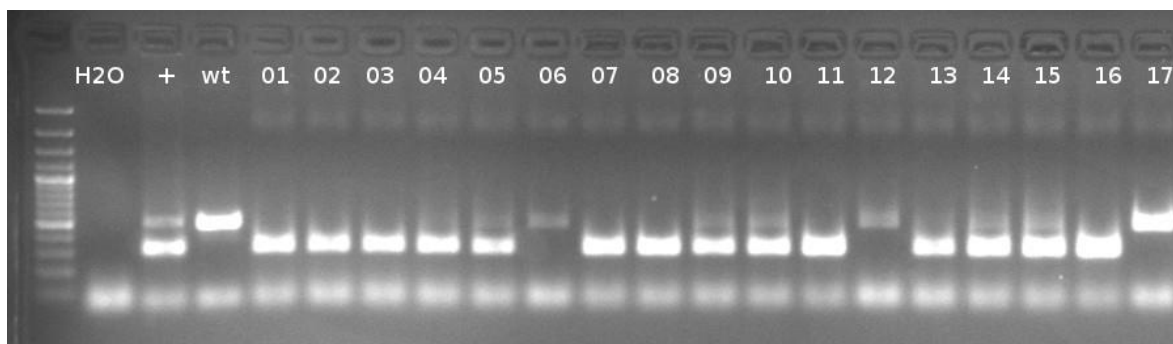

Separated by gel electrophoresis on a 2% agarose gel.