


ANK3-D1056G Genotyping Strategy

Introduction


The Idaho Technology LightScanner is a system used to perform high throughput DNA melting analysis. PCR is performed in the presence of the double stranded DNA binding dye LCGreen. After PCR samples are then heated on the LightScanner and the fluorescence emitted by bound LCGreen is monitored. As the DNA melts the LCGreen is released and so the fluorescence decreases until all the DNA has melted and all LCGreen is unbound. There are several different genotyping methods that can be used on the LightScanner.




Unlabelled probe genotyping is used to distinguish between different homozygote samples at a given SNP where scanning analysis may not have enough sensitivity. Here a 3' blocked oligonucleotide (lunaProbe) is designed that sits directly over the SNP. Asymmetric exhaustive PCR is performed using 5 times the amount of probe and opposite primer. This creates two products, one is the full PCR product between the normal primers and the other is the probe that is bound to the opposite strand.

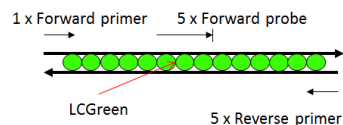


- 3' blocked oligo spanning the SNP (LunaProbe) added to PCR reaction.

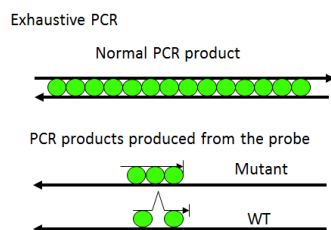

- Exhaustive, asymmetric PCR performed with unequal primer concentrations.
- Probe designed to mutant allele. Mutant allele produces the greatest fluorescence.

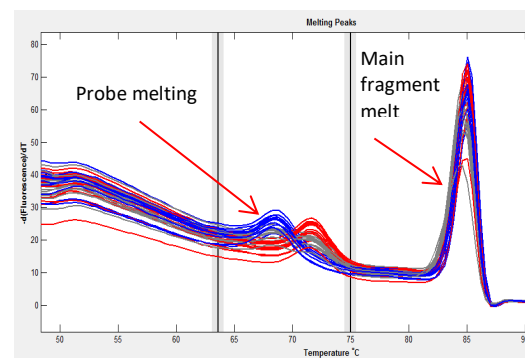

When the products are melted the probe melts at a lower temperature and by focussing analysis on this section hom, het and wt samples can be resolved.

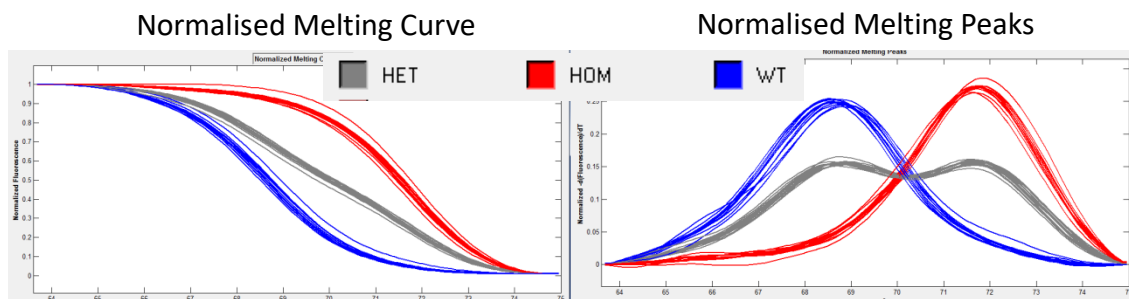



Group: Neurobehavioural
Mutation type: ENU mutant
Mutant allele: G
WT allele: A
Assay Type: LightScanner LunaProbe

Probe direction: Rev

Fragment sequence

gatacaaaaagaatcatgtggcttttttgtttgtttgtttgttttctag
ATTGTGGAAATCCCTCATTTTGGGTCCATGAGGGGGAAGGAGA GAGAACTTATCGTCCTTCG
GAGCGAGAACGGAGAGACCTGGAAG GAACATCAGTTTG G CAGTAAAAACGAAGAC CTCGCGG
AGCTT CTCAATGGCATGGATGAAG gtatcgatcccctcggatctgtggccaggacagggtgaagtgggtgaggtac

Primers/Probe sets 5'>3'

ANK3 D1056G LS FOR	GAGAACTTATCGTCCTTCGG
ANK3 D1056G LS REV	CCTTCATCCATGCCATTGAG
ANK3 D1056G PROBE rev	GTCTTCGTTTTTACTGCCAACTGATGTTC

PCR mix

HotShot master mix	5µl
LCGreen	1µl
Primer1 R(20ng/µl)	0.1µl
Primer2 F (20ng/µl)	0.5µl
Primer3 Pr R (20ng/µl)	0.5µl
DNA (1/10 dil ABI)	2µl
ddH2O	0.9µl

PCR program

LSGENO60H (annealing temperature 60 °C with hybridisation step)

Control method	Calculated
Lid control mode	Off (no need for heated lid as sample is overlaid with oil)
Lid pressure	Microplate

- | | | |
|----|------------------------|---------------|
| 1) | 95°C for 2 min | |
| 2) | 95 °C for 30 sec | PCR cycle |
| 3) | 60 °C for 30 sec | |
| 4) | 72 °C for 30 sec | |
| 5) | Cycle, step 2 55 times | |
| 6) | 95 °C for 30 sec | Hybridisation |
| 7) | 25 °C for 30 sec | |
| 8) | 15 °C for 30 sec | |



Example

