



Genotyping protocol

Rspo1 KO first with conditional
potential

Rspo1^{tm2a(KOMP)Wtsi}

IR00004991 / P4991

For any question, please contact:

Institut Clinique de la Souris - ICS - Mouse Clinical Institute

1 rue Laurent Fries, BP 10142

67404 Illkirch Cedex, France

Email: mutagenesis@igbmc.fr

Web site: <http://www-mci.u-strasbg.fr/>

TABLE OF CONTENTS

Table of contents 2

1. Genotyping protocol and data 2

 1.1. Genotyping strategy 2

 1.2. PCR protocol 4

 1.3. Picture of genotyping with various alleles 5


2. Cre and Flp genotyping method 6

1. Genotyping protocol and data

This section describes the condition used at the Mouse Clinical Institute (ICS) to genotype the **Rspo1** mouse line.

1.1. Genotyping strategy

The map below describes the position of the primers used for genotyping for each possible allele.

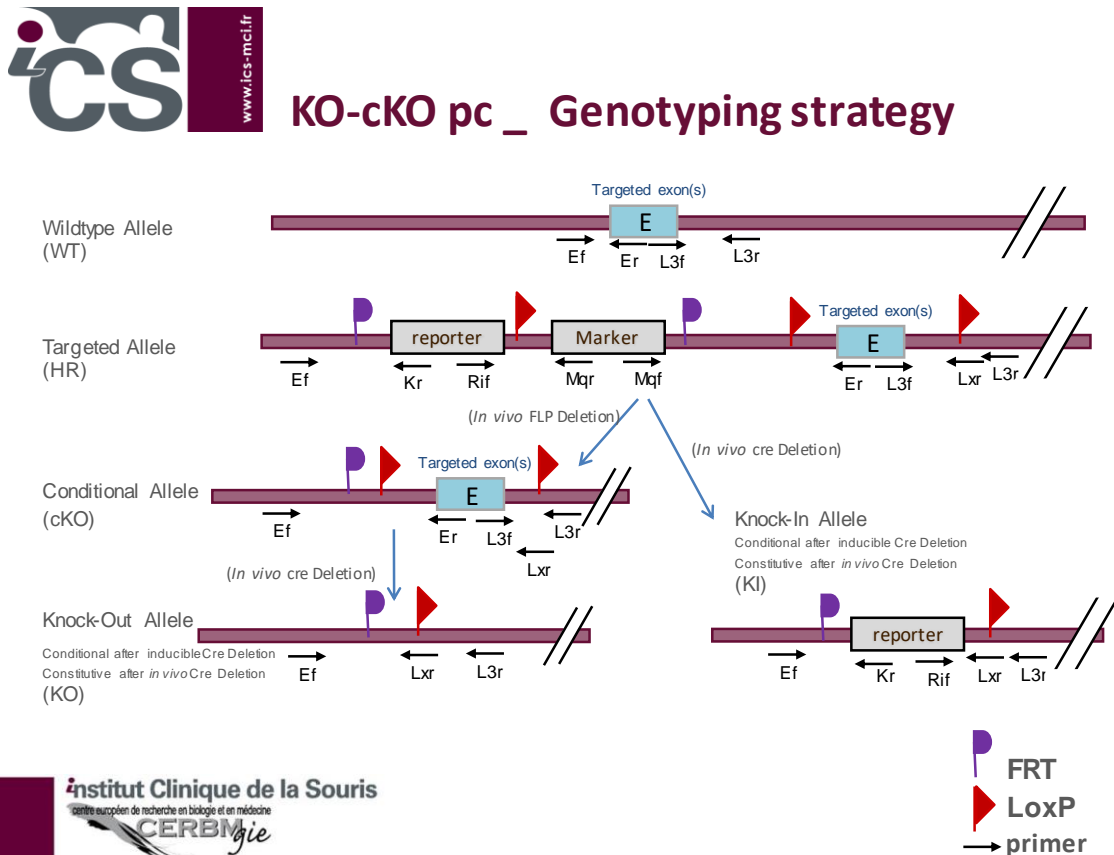

Sequence of primers used for genotyping:

Position	Sequence
Ef	GCTGTGTCAGAGTTTCTGTGGCTGG
Ef ²	CAGAGGGCCATGCACAGTCTGC
Er	GCTATCCTGGTCCTGACTAATATTGGC
Kr	CCAACAGCTTCCCCACAACGG
L3f	CTGTGATGGCTGTTGCATTTGTCCC
L3r	CTGAGCCTGGGATTGGGGAATAG
Lxr	ACTGATGGCGAGCTCAGACCATAAC
Rif	GCACATGGCTGAATATCGACGGT

²: for a selected position, a second primer was designed

PCR fragments expected size (bp):

Region analyzed	Position on the primer (see the map above)	Targeted allele (HR)	conditional allele (cKO)	KI allele	WildType allele
5' part of the selection marker	Ef / Kr	370	---	370	---
Presence of the distal loxP	L3f / L3r	278	278	---	260
Distal loxP specific PCR	L3f / Lxr	168	168	---	---
Excision of the selection marker	Ef ² / Er	7291*	387	---	243
Cre total excision	Ri1f / Lxr	3361*	---	471**	---

*: this PCR product will not be observed using our PCR genotyping conditions (see description below)

** : this PCR is only verified if mice are generated

---: no Amplicon should be obtained

1.2. PCR protocol

This section describes the composition of the mix and cycling conditions used for genotyping.

Reagents:	Volume:
- FastStart PCR Master (Roche)	7.5µl
- DNA (50ng/µl)	1.5µl
- 5' primer (100 µM)	0.06µl
- 3' primer (100 µM)	0.06µl
- Sterile H ₂ O	up to 15 µl

Cycling conditions:


Temp	Time	#Cycles
95°C	4min	1
94°C	30s	34
62°C	30s	
72°C	1min	
72°C	7min	1
20°C	5min	1

NB: These PCR conditions have been optimized for high-throughput genotyping. Adaptation to small-scale may be required.

1.3. Picture of genotyping with various alleles

Analysis of PCR products pattern was done by gel electrophoresis 2% agarose (SB buffer).

Representative genotyping picture

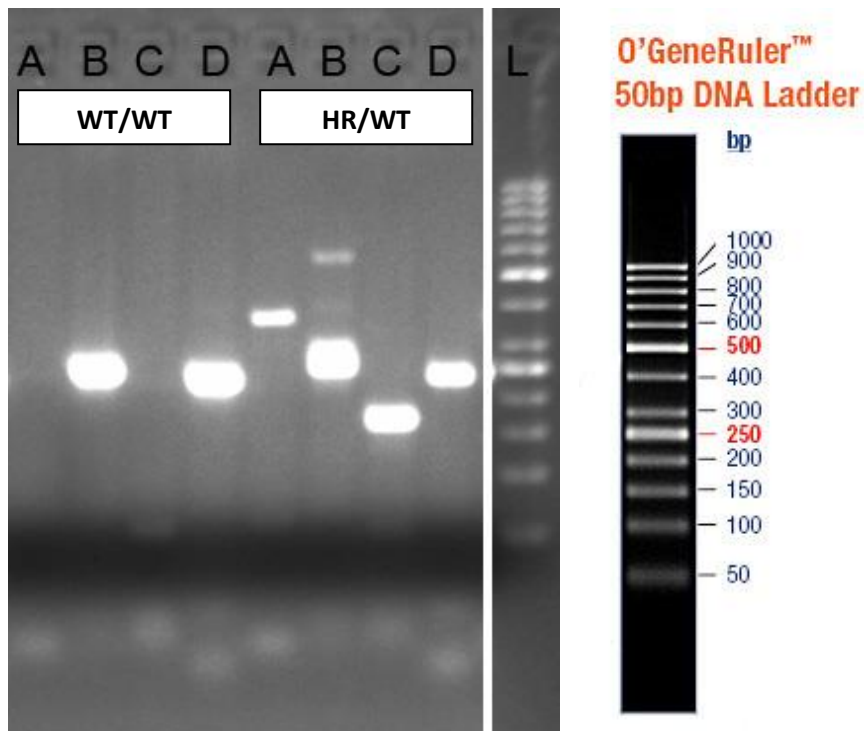

- A: 5' part of the selection marker
- B: Presence of the distal loxP
- C: Distal loxP specific PCR
- D: Excision of the selection marker
- L: O'GeneRuler 50bp DNA Ladder

2. Cre and Flp genotyping method

You will find the genotyping protocol in the publication:

[Highly-efficient, fluorescent, locus directed cre and FlpO deleter mice on a pure C57BL/6N genetic background.](#)

Birling MC, Dierich A, Jacquot S, Hérault Y, Pavlovic G.
Genesis. 2012 Jun;50(6):482-9. doi: 10.1002/dvg.20826. Epub 2012 Mar 20.